BRAUN

Age Precision™ Technology

Digital thermometer

PRT2000

GB	English	4
DE	Deutsch	12
DK	Dansk	21
ES	Español	29
FI	Suomi	38
FR	Français	46
IT	Italiano	54
NL	Nederlands	62
NO	Norsk	70
PL	Polski	78
RO	Română	86
SA	اللغة العربية	101
SE	Svenska	102
TR	Türkçe	110

Certain trademarks used under license from The Procter & Gamble Company or its affiliates.

Age Precision™ is a trademark of Helen of Troy Limited and/or its affiliates.

Kaz Europe Sàrl
A Helen of Troy Company
Q-Center, Chaux 4
CH-1030 Bussigny - Switzerland

© 2023, All rights reserved.


www.BraunHealthcare.com/uk_en/
www.BraunHealthcare.com/us_en/patents

Made and printed
in China

MD REF PRT2000EU

P/N: A008464R0
25JAN23

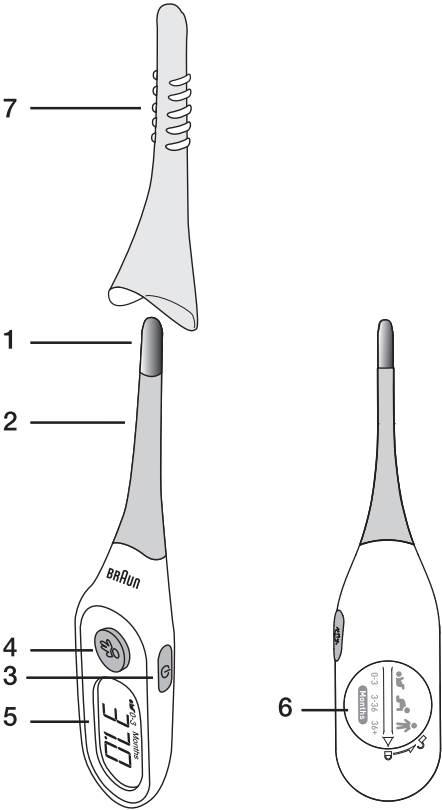Product description (See page 3)

- ① Digital sensor
- ② Comfort flex tip
- ③ Power button
- ④ Age Precision™ Technology button
- ⑤ Color-coded large display
- ⑥ Battery door (1x 3V battery)
- ⑦ Protective cap

Intended use and contraindications

The Braun thermometer is intended to be used by home users for intermittent measurement of the body temperature of individuals at the following measuring sites (oral, rectal and underarm). For home usage only. This product is not intended to diagnose any disease, but is a useful screening tool for temperature. Use of this thermometer is not intended as a substitute for consultation with your physician. The Age Precision™ Technology is not intended for pre-term babies or small-for-gestational age babies. The Age Precision™ Technology is not intended to interpret hypothermic temperatures.

Warnings and precautions

Your Braun Digital thermometer is a 3-in-1 thermometer suitable for underarm, oral or rectal use. Whichever method is used, it is recommended, for 15 minutes prior to use, to avoid eating or drinking any liquids, exercising, taking showers or baths, or smoking.

Please consult your physician if the thermometer shows elevated temperature. High, prolonged fever requires medical attention, especially for young children.

Parents/guardians should call the physician upon noticing any untoward sign(s) or symptom(s). For example, a child who exhibits irritability, vomiting, diarrhea, dehydration, changes in appetite or activity, even in the absence of fever, or who exhibits a low temperature, may still need to receive medical attention. Children who are on antibiotics, analgesics, or antipyretics should not be assessed solely on temperature readings to determine the severity of their illness.

Temperature elevation as indicated by Age Precision™ Technology may signal a serious illness, especially in neonates and infants, or adults who are old, frail, have a weakened immune system. Please seek professional advice immediately when there is a temperature elevation and if you are taking temperature for:

- Neonates and infants under 3 months. Consult your physician immediately if the temperature exceeds 37.4 °C (or 99.4 °F).
- Individuals over 60 years of age. Fever may be blunted or even absent in elderly individuals.
- Individuals having diabetes mellitus or a weakened immune system (e.g. HIV positive, cancer chemotherapy, chronic steroid treatment, splenectomy).
- Individuals who are bedridden (e.g. nursing home individual, stroke, chronic illness, recovering from surgery).
- A transplant individual (e.g. liver, heart, lung, kidney).

Do not allow children under 12 to take their temperatures unattended.

Do not allow children to walk or run during temperature taking.

This thermometer contains small parts that can be swallowed or produce a choking hazard to children. Always keep the thermometer out of children's reach.

Do not modify this equipment without authorization of the manufacturer.

Age Precision™ Technology

Clinical research shows the definition of fever changes as newborns grow into children and children grow into adults.¹ Select the appropriate age setting, take a temperature, and the display will glow green, yellow, or red to indicate whether the temperature is in a normal, fever, or high fever range, respectively. The device is only intended for the measurement of human body temperature. It is not intended to diagnose or monitor a serious medical condition or disease.

Temperature readings vary depending on the site of measurement. An oral temperature is generally 0.1 °C higher than underarm (axillary) temperature and rectal temperature is generally 0.6 °C higher than underarm temperature^{1,2}. **In Braun PRT2000 Digital thermometer, Age Precision™ Technology is optimized for underarm temperature taking.**

Temperature readings vary from person to person. The best method to determine your own normal temperature is to use the thermometer when you are feeling well. Record your reading twice a day (early morning and late afternoon). Take the average of the two temperatures. This is considered your normal body temperature. Any variation from it may indicate some sort of illness and you should consult your physician.

It is important to remember fever is a symptom, not a disease, and is just one of the many possible signs of illness. If your child exhibits signs of a serious illness, contact your physician.

How to use your Braun Age Precision™ Technology thermometer

1. Push the power button

You will hear a short beep and display segment check will be performed. After the segment check, last temperature reading will be displayed for 2 seconds. Subsequently, the age icons and age range will cycle, until the Age Precision™ Technology button is pressed.


2. Select age with Age Precision™ Technology button

Press the button to toggle through each age group.

Note: You must select an age group in order to take a measurement. Keep pressing the button until the desired age setting is displayed. The thermometer has 3 age settings:

- 0-3 months 
- 3-36 months 
- 36+ months 

The thermometer is ready to take a temperature when the display show three dashed lines.

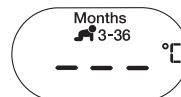3. Take temperature

Once temperature rise is detected by thermometer in measurement mode, it will begin to measure.

If the thermometer doesn't measure a temperature rise, the age icon and three dashed lines will be shown until the thermometer automatically turns off.

After measurement (generally 8 seconds after temperature change is detected), the confirmation beep indicates that an accurate temperature measurement has been taken: 1 long beep for normal temperature, 10 short beeps for elevated and high temperature. The result will be shown on the display. Interpret the reading with the help of the colour-coded display that uses medical guidelines to interpret temperature based on individual's age.

Note: Users have to shut down the thermometer and start up again to take another measurement and/or change age setting. The Braun Age Precision™ Technology thermometer turns off automatically 20 seconds after the measurement or after 90 seconds without measurement.

4. Read temperature

Braun Age Precision™ Technology colour-coded display offers a colour indication along with the temperature readout and helps you better understand the temperature reading based on patient's age. The associated colour will display 5 sec after the temperature reading was completed.

- Green backlight indicates **Normal** temperature
- Yellow backlight indicates **Elevated** temperature
- Red backlight indicates **High** temperature

In the PRT2000 digital thermometer, Age Precision™ Technology has been designed to interpret fever only for the underarm/axillary mode, according to the table below.

Site: Underarm (axillary)

Age range	Green Normal temperature	Yellow Fever	Red High fever
 0-3 Months	35.7 – 37.3 °C		> 37.3 °C
 3-36 Months	35.3 – 37.4 °C	37.5 – 38.4 °C	> 38.4 °C
 36 Months - adult	35.3 – 37.4 °C	37.5 – 39.4 °C	> 39.4 °C

If you use the thermometer in oral or rectal mode, do not take into account the color-coded display to interpret fever. Instead, refer to the oral and rectal fever thresholds, as shown in the tables below, to determine if you have fever.

Site: Oral

Age range	Green Normal temperature	Yellow Fever	Red High fever
 0-3 Months	35.8 – 37.4 °C		> 37.4 °C
 3-36 Months	35.4 – 37.6 °C	37.7 – 38.5 °C	> 38.5 °C
 36 Months - adult	35.4 – 37.7 °C	37.8 – 39.4 °C	> 39.4 °C

Site: Rectal

Age range	Green Normal temperature	Yellow Fever	Red High fever
 0-3 Months	36.3 – 37.9 °C		> 37.9 °C
 3-36 Months	35.9 – 38.1 °C	38.2 – 39.0 °C	> 39.0 °C
 36 Months - adult	35.9 – 38.2 °C	38.3 – 39.9 °C	> 39.9 °C

Note: Due to human physiology, temperature measurements taken from different sites on the same person show variability. For example, rectal readings tend to run higher than oral as well as underarm readings taken from the same person.

Temperature taking hints

Underarm/axillary mode

This method is used for babies and young children.

1. Wipe the underarm with a dry towel.
2. Place the probe tip under the arm so the tip is touching the skin and position the individual's arm next to the individual's body (Fig. 1).
With a young child, it is sometimes helpful to hug the child to keep their arm next to their body.
This ensures that the room air does not affect the reading.
3. When the peak temperature is reached, beep signals will sound. The temperature is now confirmed.
When thermometer is removed from the measuring site, the associated colour will display.

4. Turn off the thermometer and clean as recommended below.

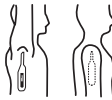
Fig. 1

Fig. 2

Oral mode

This method is suitable for children who are old enough to safely keep the thermometer inside the mouth.

1. It is important to place the tip of the thermometer well under the tongue on either side to get a good, accurate reading. Place sensor tip in areas marked with “/” under the tongue (Fig. 2). Keep patient's mouth closed and make him/her sit still to help ensure an accurate measurement.
2. When the peak temperature is reached, beep signals will sound. The temperature is now confirmed.

When thermometer is removed from the measuring site, the associated colour will display.

3. Turn off the thermometer and clean as recommended below.

Rectal mode

Commonly used for babies, young children, or when it is difficult to take an oral or underarm temperature.

1. Lubricate the thermometer's tip with a water-soluble jelly for easier insertion.
Do not use a petroleum jelly. Read the label of the lubricant to be sure you are using the right kind.
2. Lay the patient on his/her side. If the patient is an infant, the proper position for a baby is to lay on its stomach with legs hanging down, either across your knees or at the edge of a bed or changing table. This positions the infant's rectum for safe and easy insertion of the thermometer.
3. Gently insert the tip of the thermometer NO MORE THAN 1.3 cm into the rectum. If you detect resistance, stop. Hold the thermometer in place during measurement.
4. When the peak temperature is reached, beep signals will sound. The temperature is now confirmed.

When thermometer is removed from the measuring site, the associated colour will display.

5. Turn off the thermometer and clean as recommended below.

Last reading recall


When the thermometer is turned on, the display will show a test, after which the last measured temperature will be shown for approximately 1 second. Letter "M" will be displayed next to the temperature reading.

Calibration

This product has been calibrated at the time of manufacture. If used according to the use instructions, periodic calibration is not required. This device is not intended as a substitute for regular check-ups by your doctor, please consult your doctor if you have any doubt about the temperature reading.

Errors and troubleshooting

When measured temperature is higher than 42.9 °C, LCD will display "Err °C" and relevant age icon. Red backlight stays on for 5 seconds.

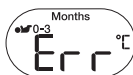When ambient temperature is out of range (see Product specifications below) or other system errors occur such as short or broken circuit, LCD will display "Err" to indicate system error. Red backlight on for 5 seconds.

Cleaning and Disinfecting


After each use, wash with soap and warm water until all visible soil is removed and let air dry. For disinfecting, wet a paper towel with >70% isopropyl alcohol, hold the tip with the wet paper towel for 60 sec, then wipe towards the tip for 3min and 30 sec. Rinse the tip with water for 30 sec and dry. Caution: Do not wet the display.

Replacing the batteries

The Braun Age Precision™ Digital thermometer comes with 1xCR1632 lithium coin battery.

When the battery symbol on LCD display twinkles continually , please replace with a new battery. The battery door design is intended to prevent a small child from easily opening the battery door to address a possible choke hazard. If the battery compartment does not close securely, stop using the thermometer and keep it away from children.

1. Insert coin into slot. Rotate door until slot is lined with UNLOCK position.
2. Replace battery.
3. Rotate battery door to LOCKED position.

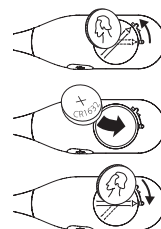BATTERY WARNING: KEEP OUT OF REACH OF CHILDREN

Swallowing can lead to chemical burns, perforation of soft tissue, and death. Severe burns can occur within 2 hours of ingestion. If you think batteries might have been swallowed or placed inside any part of the body, seek immediate medical attention immediately.

Only discard empty batteries. They should not be disposed of in the household waste, but at appropriate collection sites or at your retailer.

2 Year Limited guarantee

Please read all instructions before attempting to use this device. Please retain the receipt as proof of and date of purchase. The receipt must be presented when making any claim within the relevant guarantee period. Any claim under the guarantee will not be valid without a proof of purchase.

Your device is guaranteed for Two Years (2 Years) from date of purchase.

This guarantee covers defects in materials or workmanship that occur under normal use; defective devices meeting these criteria will be replaced free of charge.

The guarantee DOES NOT cover defects or damage resulting from abuse or failure to follow the user instructions. The guarantee becomes void if the device is opened, tampered with, or used with non-Braun branded parts or accessories, or if repairs are undertaken by unauthorized persons.

Accessories and consumables are excluded from any warranty.

For support requests, please visit www.BraunHealthcare.com/uk_en/ or find service contact information at the end of this owner's manual.

This Guarantee is applicable to Europe, Russia, Middle East, and Africa only.

UK Only: This does not affect your consumer statutory rights

Inquiry

In Europe, any serious incident (e.g., death, life-threatening injury, surgical intervention, etc.) that has occurred in relation to this thermometer should be reported to Kaz Europe (see Helen of Troy consumer service numbers) and to the competent authority of your member state.

REF, LOT, and Serial numbers of your product must be provided as they are essential to record and follow up your inquiry or claim. The LOT and Serial numbers are located in the battery compartment. The production date is given by the LOT number and can be deciphered as explained below:

The first 3 numeric digits in the LOT number represent the day of the year of manufacture. The next 2 numeric digits represent the last two numbers of the calendar year of manufacture and the letter(s) at the end designate the manufacturer of the product. (e.g.: LOT No.: 12313tav this product was made on the day 123, year 2013 at the manufacturer code tav.)

Product specifications

Displayed temperature range:	32.0 °C – 42.9 °C (89.6 °F – 109.2 °F)
Display resolution:	0.1 °C or °F
Accuracy for displayed temperature range	Maximum Laboratory Error ± 0.1 °C (± 0.2 °F) at an ambient of 23±5 °C
35.5 °C – 42 °C (95.9 °F – 107.6 °F):	±0.2 °C/0.4 °F (within 32.0 °C – 42.9 °C (89.6 °F – 109.2 °F))
Outside this range:	3V DC, 1X CR1632 lithium coin
Battery type:	500 times measurements
Battery life:	within 2 years
Service life:	2000 measurements
Operating ambient temperature range:	10 – 40 °C (50 – 104 °F) and relative humidity range of 15%–95%
Storage/Transport temperature and relative humidity:	-25 – 60 °C ± 2 °C and 15-95 RH% ± 5%
If the device is not used within specified temperature and humidity ranges the technical accuracy of the measurement cannot be guaranteed.	
This thermometer is an adjusted mode thermometer that displays the predicted temperature of the measuring site (oral or axillary or rectal) at the moment of the measure.	

Equipment with type BF applied parts

See instruction for use

Operating temperature

Storage temperature

Operating and storage humidity

Caution

Batch code

Serial number

Catalogue number

Medical device

Manufacturer

Date of manufacture

Keep it dry

Authorized representative in the European Community

Internally powered equipment.

Continuous operation.

Not intended to be sterilised.

IP22 Protected against solid foreign objects of 12.5 mm diameter and greater. Protected against vertically falling water drops when the device is tilted up to 15 degrees.

This medical device carries the CE mark and is manufactured in conformity with RoHS Directive 2011/65/EU and other applicable directive and/or regulations as specified in the EU Declaration of Conformity.

MEDICAL ELECTRICAL EQUIPMENT needs special precautions regarding EMC. For detailed description of EMC requirements please visit website (Support / Document center) or contact service center.

Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT.

Please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at your local retailer or at appropriate collection points provided in your country.

¹ Herzog L., Phillips S. Addressing Concerns About Fever. Clinical Pediatrics. 2011; 50 (#5): 383-390.

² Sund-Levander M, Forsberg C, Wahren LK. Normal oral, rectal, tympanic and axillary body temperature in adult men and women: a systematic literature review. Scand J Caring Sci 2002 June;16(2):122-8.

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